

PRODUCT DATASHEET

NOVA R5-4K

FASTCAM series by Photron



FASTCAM NOVA R5-4K

9.4-megapixel CMOS Image Sensor:

4096 x 2304 pixels (4K) at 1,250fps (NOVA R5-4K)

3840 x 2160 pixels (UHD) at 1,440fps

Maximum Frame Rate:

200,000fps (NOVA R5-4K)

Class Leading Light Sensitivity:

ISO 3,200 monochrome

ISO 640 color

Global Electronic Shutter:

2μs independent of frame rate

Dynamic Range (ADC):

12-bit monochrome

36-bit color

Compact and Lightweight:

120mm (H) x 120mm (W) x 223.2mm (D)

4.72□(H) x 4.72□(W) x 8.78□(D)

Weight: 3.5kg (7.7 lbs.)

Internal Recording Memory:

16GB, 32GB, 64GB, 128GB

Optional FASTDrive Removable High Capacity Data Storage:

4TB High-speed Solid State Drive

Fast 10-Gigabit Ethernet Interface:

Provides camera control and high-speed image download to standard PC

Fan Stop Function:

Remotely switch off cooling fans to eliminate vibration when recording at high magnifications

THE WORLDS FASTEST 4K HIGH SPEED CAMERA

The FASTCAM NOVA R5-4K high speed camera unites extraordinary CMOS image sensor innovations, and all-encompassing high-speed digital imaging expertise, to provide an ultra-HD camera with the flexibility to be used in a wide variety of applications.

True 12-bit FASTCAM NOVA R5-4K high speed camera delivers the world's fastest 4K recording rates up to 1,250 fps at 4K image resolution, and 1,440 fps at UHD image resolution. With lightning-fast shutter speeds to 2μs, blurring becomes unnoticeable, and the action appears with crisp detail. This 4k camera power is packed in a rugged, compact, lightweight camera that provides the best light-sensitivity and superior image quality.

Standard features of FASTCAM NOVA R5-4K high speed camera incorporates; an internal mechanical shutter to allow remote system calibration, a high-performance 10-Gigabit Ethernet interface for camera control and high-speed image downloading, memory segmentation that allows recording into one memory partition while downloading from another, and compatibility with a number of industry standard lens formats to allow the use of Nikon G-Type, C-mount, and Canon EF lenses. This 4K high speed camera features a sealed body design that prevents dust and corrosive particles from contaminating sensitive electronics.

The FASTCAM NOVA R5-4K high speed camera also offers an optional FASTDrive SSD which can be used for the download of images faster than 1GB per second. The intuitive and feature rich Photron FASTCAM Viewer (PFV) software is also included with each FASTCAM NOVA R5-4K camera system. Photron Device Control (SDK) allows integration of the camera with user-specific software features, and camera controlling within MATLAB® or LabView environments.

For applications that may require a high degree of motion analysis; such as automotive safety testing, weapons performance testing, fluid dynamics, and solid mechanics, will benefit from the combination of very high 4k resolution combined with high framing rates. The Photron R5-4K high speed camera falls perfectly to micro analysis applications such as Particle Image Velocimetry (PIV) or Digital Image Correlation (DIC).



Light Sensitivity:

FASTCAM NOVA R5-4K	
Monochrome models	ISO 3,200
Color models	ISO 640

Image Sensor:
FASTCAM NOVA R5-4K cameras use an advanced CMOS image sensor optimized for light sensitivity and high image quality that is unique to Photron.


A 6.5-micron pixel pitch gives a sensor size at full image resolution of 26.62 x 14.98mm (diagonal 30.55mm).

Lenses designed for both FX (35mm full frame) and also Nikon DX (digital SLR) formats are compatible with the FASTCAM NOVA R3-4K / R5-4K at UHD-4K resolution.

Sensor Type	Proprietary Design Advanced CMOS
Maximum Resolution (pixels)	4096 x 2304 pixels
Sensor Size / Diagonal	26.62 x 14.98mm / 30.55mm
Pixel Size (microns)	6.5µm x 6.5µm
Quantum Efficiency	TBD
Fill Factor	TBD
Color Matrix	Bayer CFA (single sensor)
Light Sensitivity	ISO 3,200 monochrome ISO 640 color

Shutter
Global Electronic Shutter 2.0µs independent of frame rate

FASTCAM Nova R5 Relative Spectral Response Curves – Monochrome and Color

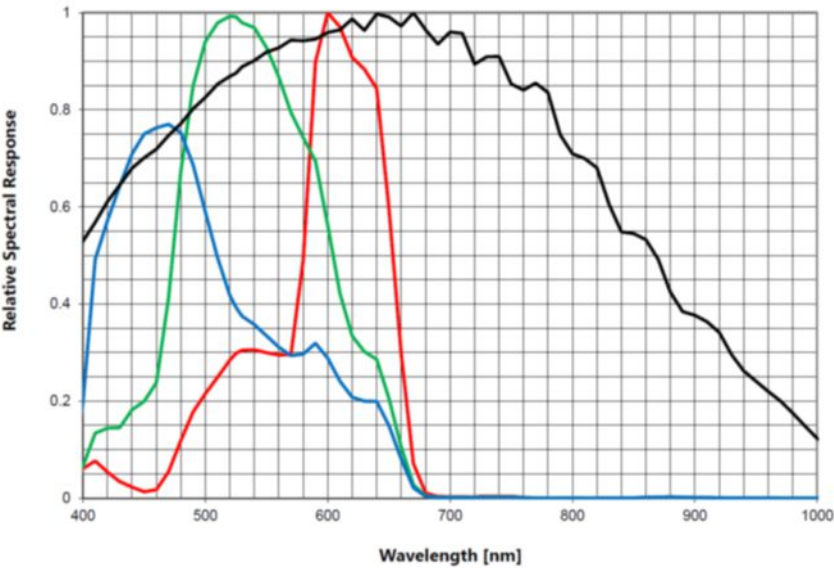

Camera Performance Specifications

Model	FASTCAM NOVA R5-4K
Full Frame Performance	1,250fps 4096 x 2304 pixels
Maximum Frame Rate	200,000fps (2048 x 8 pixels)
Minimum Exposure Time	Global electronic shutter to 2.0µs selectable independent of frame rate
Ruggedized Mechanical Calibration Shutter	Standard feature
Dynamic Range (ADC)	12-bit monochrome 36-bit color
Memory Capacity Options	16GB, 34GB, 64GB, 128GB
Memory Partitions	Up to 128 memory segments
Region of Interest	Selectable in steps of 128 pixels (horizontal) x 8 pixels (vertical)
Trigger Inputs	Selectable +/- TTL 5V and switch input (may be configured NO or NC)
Trigger Delay	Programmable on selected input / output triggers: 100ns resolution
Input / Output	Input: Trigger (TTL/Switch), sync, ready, event, IRIG Output: trigger, sync, ready, rec, exposure
Trigger Modes	Start, end, center, manual, random, random center, random manual
Time Code Input	IRIG-B (selectable at beginning or end of frame exposure)
External Sync	+/- TTL 5Vp-p Variable frequency sync
Camera Control Interface	High-speed 1/10 Gigabit Ethernet
Image Data Display	Frame rate, shutter speed, trigger mode, date/time, status, real time / IRIG time, frame count, resolution
Saved Image Formats	BMP, TIFF, JPEG, PNG, RAW, MRAW, AVI, MOV
Supported OS	Microsoft Windows operating system including: 8.1, 10 (32/64-bit)

Optional Removable Data Storage:

FASTCAM NOVA R5-4K cameras can be supplied with the Photron FASTDrive high capacity removable SSD. The ultra-high data rate FASTDrive allows a 128GB camera recording to be transferred to a removable SSD drive in approximately 2 minutes. Recorded data can then be directly accessed while coupled to the camera or the drive may be removed and inserted into the portable FASTDock station connected to any Windows PC.

High-Speed Gigabit Ethernet Interface:

FASTCAM NOVA R5-4K cameras are equipped with a high-speed 10-Gigabit Ethernet Interface to provide reliable camera control and fast download of image data.

Dedicated I/O:

A dedicated BNC connection for a contact closure hardware trigger input supporting NO, NC operation is provided. In addition, two programmable inputs and two programmable output channels provide direct connection for common tasks such as synchronization of multiple cameras and operation in conjunction with Data Acquisition (DAQ) hardware.

Ruggedized Mechanical Calibration Shutter:

The ruggedized mechanical shutter is fitted as standard to FASTCAM NOVA R5-4K cameras allows sensor black balance calibration to be carried out remotely from the system control software.

Optional Canon EF Lens Mount:

In addition to the standard C-mount and Nikon G type lens adapters, all FASTAM NOVA models support an optional Canon EF lens adapter which, through Photron FASTCAM Viewer (PFV), not only enables remote operation of lens focus and aperture but also adds Auto-Focus capability.




Camera Operation Features

Frame Synchronization	Accurate frame synchronization with other cameras and with external and unstable frequencies.
Memory Partitions	Up to 128 memory segments
Low Light Mode	Operation at minimum frame rate with separately adjustable shutter time to allow easy camera set-up and focus in ambient lighting.
Video Output	1080p live and playback via HD-SDI output
IRIG Phase Lock	Enables multiple cameras to be synchronized together with other instrumentation equipment or to a master external time source.
Internal Time Delay Generator	Allows programmable delays to be set on input and output triggers; 100ns resolution.
Event Markers	Up to ten user-entered event markers to define specific events within the recorded image sequence .
Download While Recording	FASTCAM NOVA R5-4K cameras support Partition Recording Mode, allowing image data captured in one memory partition to be downloaded while at the same time recording into another partition.
Automatic Download	The system can be set to automatically download image data to the control PC and, when download is complete to re-arm in readiness for the next trigger with automatically incremented file names.
Software Binning	Virtual pixel binning (2x2, 4x4 etc.) allows increased light sensitivity with reduced image resolution without changing camera field of view.
FASTDrive	4TB solid state drive (SSD) memory pack provides ultra high data rate transfer to removable media.

Operation Software Features

Image Calibration	2D image calibration allows the measurement of distance and angle from the image. A calibration grid overlay can be superimposed on the image.
Image Overlay	A stored reference image may be overlaid on the live image to allow accurate camera positioning to achieve the same view as a previous test.
Import of Multiple Image Sequences	Multiple image sequences can be loaded and simultaneously replayed. Timing of image sequences can be adjusted to create a common time reference. Time based synchronization allows images captured at different frame rates to be synchronized.
High Dynamic Range Mode	Making use of the full sensor dynamic range, HDR mode allows enhanced detail in both light and dark areas of an image to be displayed simultaneously.
Background Subtraction	In order to highlight subtle changes in an image, Background Subtraction allows a reference image to be subtracted from a recorded sequence. Details including propagation of shock waves and surface changes during impact can be visualized using the feature.
Line Profile	A line profile representing grey levels along a line drawn across any region of the image is displayed. In live mode the Line Profile can be used to ensure optimum image focus is achieved.
Histogram	A histogram displaying grey levels within a user-defined image area is displayed. In live mode the Histogram can be used to ensure that optimum exposure levels are set for the scene being recorded.

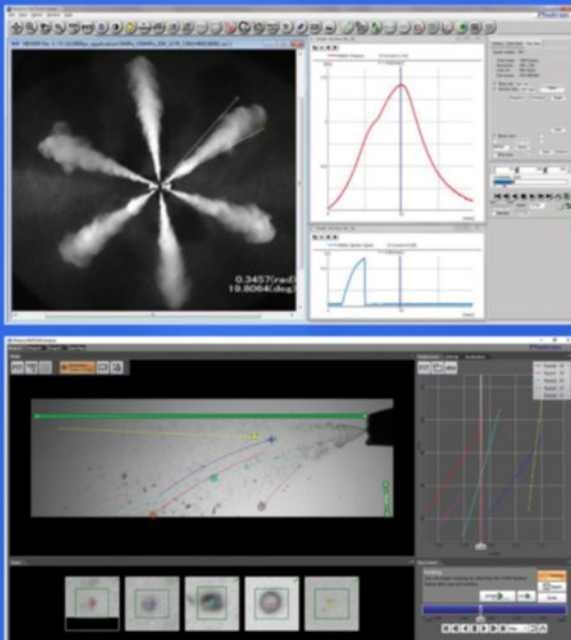

Photron FASTCAM Viewer:

Photron FASTCAM Viewer software (PFV) has been designed to provide an intuitive and feature rich user interface for the control of Photron high-speed cameras, data saving, image replay and simple motion analysis. Advanced operation menus provide access to features for advanced camera operation and image enhancement. Tools are provided to allow image calibration and easy measurement of angles and distances from image data. Also included are a C++ SDK and wrappers for LabView and MATLAB ®.

An optional software plug-in module provides synchronization between Photron high-speed cameras and data acquired through National Instruments data acquisition systems. Synchronized data captured by the DAQ system provides waveform information which can be viewed alongside high-speed camera images.

Photron FASTCAM Analysis:

PFV software allows image sequences to be exported directly to optional Photron FASTCAM Analysis (PFA) Motion Analysis software. This entry level Motion Analysis software with an on screen 'step by step guide' function provides automated tracking of up to 5 points using feature or correlation tracking algorithms for the automated analysis of motion within an image sequence.

FASTCAM NOVA R5-4K	
Resolution	Maximum Frame Rate
(h x v pixels)	NOVA R5-4K
4096 x 2304	1,250
4096 x 1504	2,000
4096 x 1024	2,880
4096 x 512	5,500
4096 x 256	11,250
4096 x 128	20,000
4096 x 64	40,000
4096 x 8	158,400
3840 x 2160	1,440
2048 x 2048	2,880
2048 x 1024	5,500
2048 x 512	11,250
2048 x 256	20,000
2048 x 128	40,000
2048 x 64	72,000
2048 x 8	200,000
1920 x 1080	5,280
1280 x 720	8,000
1024 x 1024	5,500
640 x 480	12,000
512 x 512	11,250

* Specifications subject to change without notice.

Variable Region of Interest:

Region of Interest (ROI) or sub-windowing allows a user-specified portion of the sensor to be defined to capture images. By using a reduced portion of the image area, the frame rate at which images are recorded can be increased. FASTCAM NOVA R5-4K cameras allow the ROI to be set in increments of 128 pixels horizontal and 8 pixels vertical.

External Frame Synchronization:

FASTCAM NOVA R5-4K cameras can be fully synchronized with an external source to allow the timing of when each individual image is captured to be precisely referenced. The camera can be accurately synchronized to unstable frequencies allowing complex events such as combustion in rapidly accelerating or decelerating engines to be recorded and studied.

Record During Download Operation:

FASTCAM NOVA R5-4K cameras recording memory can be divided into multiple active sections. The user can record an on-going event in one memory partition while at the same time downloading a previously recorded image sequence in order to improve workflow and optimize camera operation.



Mechanical and Environmental Specifications

Mechanical	
Lens Mount	M42, F-mount (G-type lens compatible) and C-mount provided - Optional lens mounts available include Canon EF remote control mount
Camera Mountings	3/8 - 16 UNC, 1/4 - 20 UNC & 4 x M6 (base and side), 2 x 1/4 - 20 UNC (top)
External Dimensions	
Camera Body (excluding protrusions)	120mm (H) x 120mm (W) x 223.2mm (D) 4.72" (H) x 4.72" (W) x 8.78" (D)
Weight	
Camera Body	3.5kg (7.2lbs)
Environmental	
Operating Temperature	-10 to 45C, 14° to 113°F
Storage Temperature	-20 to 60C, -4° to 140°F
Humidity	85% or less (non-condensing)
Cooling	Internal fan cooling (fan-off mode supported)
Operational Shock	30G, 11ms, 6-axes 10 times/axis
Power	
AC Power (with supplied adapter)	100 to 240V, 50 to 60Hz
DC Power (primary input)	22 to 32V, 120VA
DC Power (battery input)	22 to 32V, 120VA



Nikon G-Type Compatible Lens Mount:

FASTCAM NOVA R5-4K cameras are equipped with an objective lens mount compatible with readily available Nikon G-type lenses. Controls provided within the lens mount allow the control of lens aperture on lenses without external iris control.

Optional Canon EF Lens Mount:

An optional lens mount supporting Canon EF lenses is available for remote control of lens aperture and focus including Auto-Focus capability through Photron PFV software.

Operation Environments:

The 'sealed body' design of FASTCAM NOVA R5-4K camera ensure optimum air flow and prevents dust and corrosive particles from being ingested within the internal camera body where they can damage sensitive electronics. The fans may be disabled during recording for any vibration sensitive measurements.

FASTCAM NOVA R5-4K camera have been extensively tested to ensure operation for extended periods in ambient temperatures up to 45 degrees C.

Fan Stop Function:

Remotely switch off cooling fans to eliminate vibration when recording at high magnifications.

Specifications subject to change without notice.

PHOTRON NORTHEAST DISTRIBUTOR

Experts in Video
Instrumentation

TECH IMAGING

SERVICES, INC.

imaging@techimaging.com

978-740-0063

WWW.TECHIMAGING.COM

FOR MORE INFORMATION , PRICING, OR AVAILABILITY, PLEASE CONTACT US AT 978-740-0063